Guardian[™] Link (3)

G

Transmitter

Contacts:

Africa: Medtronic Africa (Pty) Ltd. Tel: +27 (0) 11 677 4800

Argentina: Corpomedica S.A. Tel: +(11) 4 814 1333 Medtronic Directo 24/7: +0800 333 0752

Australia: Medtronic Australasia Pty. Ltd. Tel: 1800 668 670

Azerbaijan: Isomed Tel: +994 (12) 464 11 30

Bangladesh: Sonargaon Healthcare Pvt Ltd. Mobile: (+91)-9903995417 or (+880)-1714217131

Belarus: ОДО "Баджин" Tel: +375 17 313 0990

België/Belgique: N.V. Medtronic Belgium S.A. Tel: 0800-90805

Bosnia and Herzegovina: Medimpex d.o.o. Tel: +387 33 476 444 or +387 33 476 400 Fax: +387 33 476 401 or +387 33 432 241

Brasil: Medtronic Comercial Ltda. Tel: +(11) 2182-9200 Medtronic Directo 24/7: +0800 773 9200

Bulgaria: RSR Ltd. Tel: +359 885 428 900

Canada: Medtronic of Canada Ltd. Tel: 1-800-284-4416 (toll free/sans-frais)

Chile: Medtronic Chile Tel: +(9) 66 29 7126 Medtronic Directo 24/7: +1 230 020 9750 Medtronic Directo 24/7 (From Santiago): +(2) 595 2942

China: Medtronic (Shanghai) Ltd. 24 Hour Help (Cell): +86 400-820-1981 24 Hour Help (Landline): +86 800-820-1981 Colombia: Medtronic Latin America Inc. Sucursal Colombia Tel: +(1) 742 7300 Medtronic Directo 24/7 (Landline): +01 800 710 2170 Medtronic Directo 24/7 (Cellular): +1 381 4902

Croatia: Medtronic Adriatic d.o.o. Tel: +385 1 488 11 20 Fax: +385 1 484 40 60

Danmark: Medtronic Danmark A/S Tel: +45 32 48 18 00

Deutschland: Medtronic GmbH Geschäftsbereich Diabetes Telefon: +49 2159 8149-370 Telefax: +49 2159 8149-110 24-Stdn-Hotline: 0800 6464633

Eire: Accu-Science LTD. Tel: +353 45 433000

España: Medtronic Ibérica S.A. Tel: +34 91 625 05 42 Fax: +34 91 625 03 90 24 horas: +34 900 120 330

Europe: Medtronic Europe S.A. Europe, Middle East and Africa Headquarters Tel: +41 (0) 21-802-7000

France: Medtronic France S.A.S. Tel: +33 (0) 1 55 38 17 00

Hellas: Medtronic Hellas S.A. Tel: +30 210677-9099

Hong Kong: Medtronic International Ltd. Tel: +852 2919-1300 To order supplies: +852 2919-1322 24-hour helpline: +852 2919-6441

India: India Medtronic Pvt. Ltd Tel: (+91)-80-22112245 / 32972359 Mobile: (+91)-9611633007 Patient Care Helpline: 1800 209 6777

Indonesia: Medtronic International Ltd. Tel: +65 6436 5090 or +65 6436 5000 **Israel:** Medtronic World Trade Corporation Tel: +972 9972 4400

Italia: Medtronic Italia S.p.A. Tel: +39 02 24137 261 Fax: +39 02 24138 210 Servizio assistenza tecnica: Nº verde: 800 60 11 22

Japan: Medtronic Japan Co. Ltd. Tel: +81-3-6776-0019 24 Hr. Support Line: 0120-56-32-56

Kazakhstan: Medtronic BV in Kazakhstan Tel: +7 727 311 05 80 (Almaty) Tel: +7 717 224 48 11 (Astana) Круглосуточная линия поддержки: 8 800 080 5001

Latin America: Medtronic, Inc. Tel: 1(305) 500-9328 Fax: 1(786) 709-4244

Latvija: Ravemma Ltd. Tel: +371 7273780

Macedonia: Alkaloid Kons Dooel Tel: +389 2 3204 430

Magyarország: Medtronic Hungária Kft. Tel: +36 1 889 0688

Malaysia: Medtronic International Ltd. Tel: +603 7946 9000

Middle East and North Africa: Regional Office Tel: +961-1-370 670

Montenegro: Glosarij Tel: +382 20 642 495 Fax: +382 20 642 540

México: Medtronic Servicios S. de R. L. de C.V. Tel (México DF): +(11) 029 058 Tel (Interior): +01 800 000 7867 Medtronic Directo 24/7 (from México DF): +(55) 36 869 787 Medtronic Directo 24/7: +01 800 681 1845

Nederland, Luxembourg: Medtronic B.V. Tel: +31 (0) 45-566-8291 Gratis: 0800-3422338 New Zealand: Medica Pacifica Phone: 64 9 414 0318 Free Phone: 0800 106 100

Norge: Medtronic Norge A/S Tel: +47 67 10 32 00 Fax: +47 67 10 32 10

Philippines: Medtronic International Ltd. Tel: +65 6436 5090 or +65 6436 5000

Россия: ООО «Медтроник» Tel: +7 495 580 73 77 Круглосуточная линия поддержки: 8 800 200 76 36

Polska: Medtronic Poland Sp. z o.o. Tel: +48 22 465 6934

Portugal: Medtronic Portugal Lda Tel: +351 21 7245100 Fax: +351 21 7245199

Puerto Rico: Medtronic Puerto Rico Tel: 787-753-5270

Republic of Korea: Medtronic Korea, Co., Ltd. Tel: +82.2.3404.3600

Romania: Medtronic BV Reprezentanta Tel: +40 372 188 000

Schweiz: Medtronic (Schweiz) AG Tel: +41 (0)31 868 0160 24-Stunden-Hotline: 0800 633333 Fax Allgemein: +41 (0)318680199

Serbia: Medtronic B.V. Serbia Tel: +381 11 2095 900

Singapore: Medtronic International Ltd. Tel: +65 6436 5090 or +65 6436 5000

Slovenija: Zaloker & Zaloker d.o.o. brezplačna številka: 080 1880 Tel: +386 1 542 51 11

Slovenská republika: Medtronic Slovakia, s.r.o. Tel: +421 26820 6942 HelpLine: +421 26820 6986 Sri Lanka: Swiss Biogenics Ltd. Mobile: (+91)-9003077499 or (+94)-777256760

Suomi: Medtronic Finland Oy Tel: +358 20 7281 200 Help line: +358 400 100 313

Sverige: Medtronic AB Tel: +46 8 568 585 20 Fax: +46 8 568 585 11

Taiwan: Medtronic (Taiwan) Ltd. Tel: 02-21836000 Toll free: +886-800-005285

Thailand: Medtronic (Thailand) Ltd. Tel: +662 232 7400

Türkiye: Medtronic Medikal Teknoloji Ticaret Ltd. Sirketi. Tel: +90 216 4694330

USA: Medtronic Diabetes Global Headquarters 24 Hour HelpLine: +1-800-646-4633 To order supplies: +1-800-843-6687

Ukraine: Medtronic B.V. Representative office in Ukraine Tel: +38 044 392 04 01 Лінія цілодобової підтримки: 0 800 508 300


United Kingdom: Medtronic Ltd. Tel: +44 1923-205167

Österreich: Medtronic Österreich GmbH Tel: +43 (0) 1 240 44-0 24 – Stunden – Hotline: 0820 820 190

Česká republika: Medtronic Czechia s.r.o. Tel: +420 233 059 111 Non-stop helpLine (24/7): +420 233 059 059 Zákaznický servis (8:00 - 17:00): +420 233 059 950

English

The Guardian[™] Link (3) transmitter is a component of the continuous glucose monitoring (CGM) system for the MiniMed[®] 670G insulin pump. The transmitter is compatible only with the Guardian[™] Sensor (3) glucose sensor. The transmitter collects data from the sensor. The transmitter then wirelessly sends the data to the insulin pump.

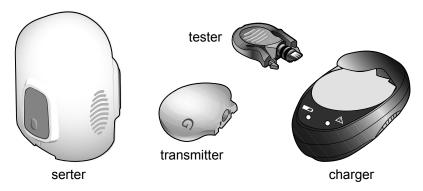

Guardian Link (3) transmitter kit components

A complete transmitter kit includes the following components:

- Guardian Link (3) transmitter (MMT-7811)
- Charger (MMT-7715)

Two testers (MMT-7736L)

• One-press serter (MMT-7512)

Indications for use

The Guardian Link (3) transmitter is intended for use with MiniMed 670G System. The transmitter powers the glucose sensor, collects and calculates sensor data, and wirelessly sends the data to the MiniMed 670G pump. The transmitter is intended for single-patient, multi-use.

Contraindications

None known.

Warnings

• Do not use the transmitter adjacent to other electrical equipment which may cause interference with the normal system operation. This includes mobile communication devices such as cell phones, GPS navigation systems, and other devices that have an output transmitter power greater than 1W. Other electrical

equipment that may compromise normal system operation has been contraindicated. For more information, see *Exposure to magnetic fields and radiation, on page 2*.

- Always refer to the sensor user guide for all precautions, warnings, and instructions relating to the sensor. Not referring to the sensor user guide can result in serious injury or damage to the sensor.
- Do not allow children to put small parts in their mouth. This product poses a choking hazard for young children.
- Do not change or modify the device unless expressly approved by Medtronic Diabetes. Modifying the device can cause serious injury, interfere with your ability to operate the device, and void your warranty.
- Do not use the tester if it comes in contact with blood. Touching blood can cause infection. Dispose of the tester according to the local regulations for medical waste disposal, or contact your healthcare professional for disposal information.
- Bleeding may occur after inserting the sensor. Always make sure that the site is not bleeding before connecting the transmitter to the sensor. Blood can get into the transmitter connector and damage the device. Discard the device if damaged. If bleeding occurs, apply steady pressure with a sterile gauze or clean cloth at the insertion site until bleeding stops. After bleeding stops, connect the transmitter to the sensor.
- Contact the 24 Hour HelpLine if you experience any adverse reactions associated with the transmitter or sensor. Adverse reactions can cause serious injury.

Exposure to magnetic fields and radiation

- Do not expose your transmitter to Magnetic Resonance Imaging (MRI) equipment, diathermy devices, or other devices that generate strong magnetic fields (for example, x-ray, CT scan, or other types of radiation). Exposure to a strong magnetic field has not been evaluated and can cause the device to malfunction, result in serious injury or be unsafe. If your transmitter is exposed to a strong magnetic field, discontinue use and contact the 24 Hour HelpLine for further assistance.
- Always remove your sensor and transmitter before entering a room that has x-ray, MRI, diathermy, or CT scan equipment. Exposure to a strong magnetic field has not been evaluated and can cause the device to malfunction, result in serious injury or be unsafe. If your sensor or transmitter is exposed to a strong magnetic field, discontinue use and contact the 24 Hour HelpLine for further assistance.
- Always carry the Medical emergency card provided with your device when you are traveling. The Medical emergency card provides critical information about airport security systems, and using your transmitter on an airplane that can help you and others. Not following the guidance on the Medical emergency card could result in serious injury.

Precautions

- Only use the Guardian Sensor (3) (MMT-7020) glucose sensor with the transmitter. Do not use any other sensor. Other sensors are not intended for use with the transmitter, and will damage the transmitter and the sensor.
- Only use the green colored tester (MMT-7736L) with the transmitter. Pockets on the transmitter are visible when connected to the tester. Do not use any other test plug. Other test plugs are not intended for use with the transmitter, and will damage the transmitter and the tester.
- Always use the tester when cleaning the transmitter. Do not use any other test plug with the transmitter. Use of another test plug can allow water to get into the transmitter or can prevent proper cleaning. Water can damage the transmitter.
- Do not twist the tester or sensor while attached to the transmitter. Twisting the tester or sensor will damage the transmitter.
- Do not allow the tester to come in contact with any liquid when not connected to the transmitter. A wet tester can damage the transmitter.
- Do not allow the transmitter to come in contact with any liquid when not connected to a sensor or to the tester. Moisture will damage the transmitter and a wet transmitter can damage the sensor.
- Do not clean the o-rings on the tester with any substances. Cleaning the o-rings can damage the tester.

Radio Frequency (RF) communication

This device complies with the United States Federal Communications Commission (FCC) and international standards for electromagnetic compatibility. This device complies with Part 15 of the FCC Rules. Operation is subject to two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

This device has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This device generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this device does cause harmful interference to radio or television reception, which can be determined by turning the device off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

• Reorient or relocate the receiving antenna.

- Increase the separation between the device and the receiver.
- Decrease the distance between the transmitter and the insulin pump to 6 feet (1.8 meters) or less.
- Increase the separation between the transmitter and the equipment that is receiving or emitting interference.
 - **Note:** Harmful interference is defined by the FCC as follows. Any emission, radiation or induction that endangers the functioning of a radio navigation service or of other safety services or seriously degrades, obstructs or repeatedly interrupts a radio communications service operating in accordance with FCC rules.

Changes or modifications made to this equipment not expressly approved by Medtronic Diabetes could void the user's authority to operate the equipment.

Directive 1999/5/EC

Medtronic declares that this product is in conformity with the essential requirements of Directive 1999/5/EC on Radio and Telecommunications Terminal Equipment.

For additional information, contact Medtronic MiniMed at the address or phone number provided on the back cover.

IEC60601-1-2:2007; Special EMC Precautions for Medical Electrical Equipment

- Special Precautions regarding Electromagnetic Compatibility (EMC): This body worn device is intended to be operated within a reasonable residential, domestic, public or work environment, where common levels of radiated "E" (V/m) or "H" fields (A/m) exist; such as cellular phones, WiFi, Bluetooth, electric can openers, microwave and induction ovens. This device generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the provided instructions, may cause harmful interference to radio communications.
- 2 Portable and mobile RF communications equipment can affect Medical Electrical Equipment as well. If you encounter RF interference from a mobile or stationary RF transmitter, move away from the RF transmitter that is causing the interference.

Assistance

Medtronic MiniMed provides a 24 Hour HelpLine for assistance. When calling the HelpLine, please have the serial number of your device available. The serial number and the 24 Hour HelpLine phone number are listed on the back of your device.

Department	Telephone number
24 Hour HelpLine (calls within the United States)	800 646 4633
24 Hour HelpLine (calls outside the United States)	+1 818 576 5555
Website	www.medtronicdiabetes.com

Charger


The transmitter contains a non-replaceable, rechargeable battery that you can recharge as needed with the charger. The charger has a green light that shows the charging status and a red light that communicates any problems during charging. If you see a red light, see the Troubleshooting section. The charger needs one AAA alkaline battery.

Note: If the battery is installed incorrectly or is low, the charger will not work. Repeat the battery installation steps using a new battery.

Installing a battery in the charger

To install a battery in the charger:

- 1 Push the battery cover in and slide it off (as shown in the image in step 3).
- 2 Insert a new AAA alkaline battery. Make sure the + and symbols on the battery align with these same symbols shown on the charger.
- 3 Slide the cover back on the charger until it clicks into place.

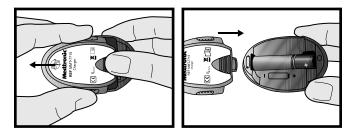

Charging the transmitter


Caution: Always charge the transmitter before inserting your sensor. A depleted transmitter does not function. A fully charged transmitter works at least seven days without recharging. A depleted transmitter can take up to two hours to recharge.

Caution: Do not store the transmitter on the charger for more than 60 days. Disconnect and reconnect to the charger to re-charge again before use. If the transmitter is left on the charger for more than 60 days, the battery will be permanently damaged.

To charge the transmitter:

- 1 Push the two components together to connect the transmitter to the charger.
- 2 Within 10 seconds after the transmitter is connected, a green light on the charger will flash for one to two seconds as the charger powers on. For the rest of the charging time, the green light on the charger will continue to flash in a pattern of four flashes with a pause between the four flashes.
- 3 When charging is complete, the green light on the charger will stay on, without flashing, for 15 to 20 seconds and then turn off.
- 4 After the green charger light turns off, disconnect the transmitter from the charger. The green light on the transmitter will flash 10 times and then turn off.



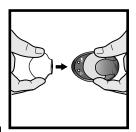

Pairing your transmitter

Always refer to the system user guide for instructions on how to pair your transmitter to your pump.

Inserting the sensor

Always refer to your sensor user guide for instructions on how to insert the sensor.

Connecting the transmitter to the sensor

Before proceeding, have your system user guide available.

To connect the transmitter to the sensor:

- 1 After the sensor is inserted, consult your sensor user guide for details on applying the required tape before connecting the transmitter.
- 2 Hold the rounded end of the inserted sensor to prevent it from moving during connection.

- 3 Hold the transmitter as shown. Line up the two notches on the transmitter with the side arms of the sensor. The flat side of the transmitter should face the skin.
- 4 Slide the transmitter onto the sensor connector until the sensor arms snap into the notches on the transmitter. If the transmitter is properly connected, and if the sensor has had enough time to become hydrated, the green light on the transmitter will flash 6 times.

Note: If the transmitter does not flash, see Troubleshooting, on page 14.

- 5 When the transmitter light flashes green after connecting to the sensor, use your pump to start the sensor. For more instructions, see your system user guide.
- 6 Attach the adhesive tab of the sensor to the transmitter.
- 7 Refer to the sensor user guide for instructions on how to apply a second tape.
- 8 Follow the instructions that appear on the pump screen or in your system user guide.

Disconnecting the transmitter from the sensor

Before proceeding, have your system user guide available.


To disconnect the transmitter from the sensor:

- 1 Carefully remove any tape from the transmitter and sensor.
- 2 Remove the adhesive tab from the top of the transmitter.
- 3 Hold the transmitter as shown, and pinch the flexible side arms of the sensor between your thumb and forefinger.
- 4 Gently pull the transmitter away from the sensor.
- 5 Follow the instructions that appear on the pump or in your system user guide.

Removing the sensor

Always refer to the sensor user guide for instructions on how to remove the sensor.



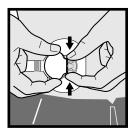

Bathing and swimming

After the transmitter and sensor are connected, they form a waterproof seal to a depth of 8 feet (2.4 meters) for up to 30 minutes. You can shower and swim without removing them.

Tester


The tester is used to test the transmitter to make sure it is working. It is also used as a required component for cleaning the transmitter. Properly connecting the tester to the transmitter ensures that fluids do not come in contact with the connector pins inside the transmitter. Fluids can cause connector pins to corrode and affect the performance of the transmitter.

Do not twist the tester while attached to the transmitter. This will damage the transmitter.

The tester can be used for one year. If you continue to use the tester for more than one year, the connector pins inside the transmitter could be damaged, because the tester cannot continue to provide a waterproof seal. For instructions on how to check the connector pins, see *Inspecting the transmitter connector pins, on page 9*.


Caution: Only use the green colored tester (MMT-7736L) with the transmitter. Pockets on the transmitter are visible when connected to the tester. Do not use any other test plug. Other test plugs are not intended for use with the transmitter, and will damage the transmitter and the tester.



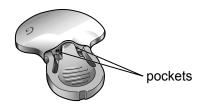

Inspecting the transmitter connector pins

This image is an example of how the connector pins should look.

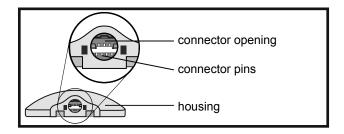

Look inside the connector opening of the transmitter to make sure that the connector pins are not damaged or corroded. If the connector pins are damaged or corroded, the transmitter cannot communicate with the charger or pump. Contact the 24 Hour HelpLine. It may be time to replace your transmitter.

Also look for moisture inside the connector opening. If you see any moisture, allow the transmitter to dry for at least one hour. Moisture inside the connector opening could cause the transmitter to not work properly, and could cause corrosion and damage over time.


Connecting the tester for testing or cleaning

Before proceeding, have your system user guide available.

To connect the tester:

- 1 Hold the transmitter and the tester as shown. Line up the flat side of the tester with the flat side of the transmitter.
- 2 Push the tester into the transmitter until the flexible side arms of the tester click into the notches on both sides of the transmitter.

When properly connected, the green light on the transmitter flashes 6 times.

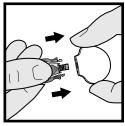

- 3 To test the transmitter, check the sensor icon on the pump to ensure that the transmitter is sending a signal (see your system user guide).
- 4 To clean the transmitter, see *Cleaning the transmitter, on page 10*.
- 5 After testing or cleaning, disconnect the tester from the transmitter.

Disconnecting the tester

To disconnect the tester:

- 1 Hold the transmitter body as shown and pinch the side arms of the tester.
- 2 With the tester arms pinched, gently pull the transmitter away from the tester.

Note: To save transmitter battery life, do NOT leave the tester connected after cleaning or testing.

Cleaning the transmitter

The transmitter is a single-patient use device and not intended for multi-patient use.

Warning: Do not discard the transmitter in a medical waste container or expose it to extreme heat. The transmitter contains a battery that may ignite, and result in serious injury.

Note: The tester is a required component for cleaning the transmitter. For details, see Tester, on page 8.

Always clean the transmitter after each use.


To clean the transmitter, you need the following materials:

- mild liquid soap (for example, Ivory[®] liquid soap)
- soft-bristled toddler toothbrush
- container
- clean, lint-free dry cloths

You can find these supplies at Walmart, Target, or http://www.amazon.com/.

Use life


The transmitter can be cleaned up to 122 times or one year, whichever comes first. Discard the transmitter at this point. If you continue to use the transmitter beyond 122 times or one year, the cleaning process may damage the device. Contact Medtronic to order a new transmitter.

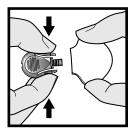


Warning: Do not use the device if you see any cracking, flaking, or damage to the housing. Cracking, flaking, or damage to the housing are signs of deterioration. Deterioration of the housing can affect the ability to properly clean the transmitter, and result in serious injury. Call the 24 Hour HelpLine and discard the device according to local regulations for battery disposal (nonincineration), or contact your healthcare professional for disposal information.

To clean the transmitter:

- 1 Wash your hands thoroughly.
- 2 Attach the tester to the transmitter.

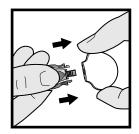


- 3 If there is adhesive residue on the transmitter, see *Removing adhesive residue, on page 13.*
- 4 Rinse the transmitter under room temperature tap water for at least one minute, and until visibly clean. Make sure all hard-to-reach areas are rinsed completely.

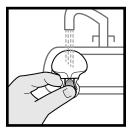

5 Prepare a mild liquid soap solution using 1 teaspoon (5 milliliters) of mild liquid soap per 1 gallon (3.8 liters) of room temperature tap water.

6 With the tester still attached, submerge the transmitter in the mild liquid soap solution and soak for one minute.

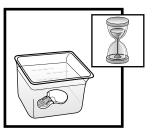

7 Holding the tester, brush the entire surface of the transmitter using a soft-bristled toddler toothbrush. Make sure to brush all hard-to-reach areas until visibly clean.




8 Rinse the transmitter under running room temperature tap water for at least one minute, and until all visible liquid soap is gone.



9 Dry the transmitter and tester with a clean, dry cloth.



- 10 Place the transmitter and tester on a clean, dry cloth and air dry them completely.
- 11 Disconnect the tester from the transmitter by gently squeezing the arms of the tester.

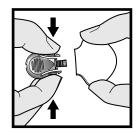

Removing adhesive residue

You may need to perform this procedure if there is adhesive residue present on the transmitter. If you visually inspect the transmitter and see adhesive residue on it, follow the instructions below.

To remove adhesive residue, you need Detachol[®] medical adhesive remover, and cotton swabs. You can buy Detachol at http://www.amazon.com/, http:// www.medtronicdiabetes.com, or by calling 800 646 4633.

Note: During testing, Medtronic MiniMed used Detachol to remove the adhesive residue from the transmitter.

To remove adhesive residue:

1 Make sure the tester is attached to the transmitter.

2 Holding the tester, saturate a cotton swab in the Detachol solution and gently rub the adhesive residue on the transmitter until it is fully removed.



3 Continue with the cleaning procedure. See *Cleaning the transmitter, on page 10* for details.

Cleaning the charger

This procedure is for general cleaning as required, based on physical appearance.

- Caution: Do not immerse the charger in water or any other cleaning agent. The charger is not waterproof. Water can damage the charger, and cause the device to malfunction.
- Warning: Dispose the charger according to the local regulations for battery disposal, or contact your healthcare professional for disposal information. The charger may ignite upon incineration.

To clean the charger:

- 1 Wash your hands thoroughly.
- 2 Use a damp cloth with mild cleaning solution, such as a dishwashing detergent, to clean any dirt or foreign material from the outside of the charger. Never use organic solvents, such as paint thinner or acetone, to clean the charger.
- 3 Place the charger on a clean, dry cloth and air dry for two to three minutes.

Troubleshooting

The following table contains troubleshooting information for the transmitter, charger, and tester. For more information about troubleshooting, see your system user guide.

Problem	Likely Cause(s)	Resolution
You connected the transmitter to the charger and no lights came on.	The transmitter connec- tor pins are damaged or corroded. Your charger battery has no power or no battery is inserted.	1 Check the transmitter connector pins for damage or corrosion. For more informa- tion about your connector pins, see <i>In-</i> <i>specting the transmitter connector pins,</i> <i>on page 9.</i> If the pins are damaged or cor- roded, contact the 24 Hour HelpLine. It may be time to replace your transmitter.
		2 If there is no damage to the connector pins, replace the battery in the charger. For instructions on replacing your charger battery, see <i>Installing a battery in the</i> <i>charger, on page 5.</i>
During charging, the flashing green light on the charger turns off and you see a longer flashing red light on the charger.	Your charger battery is low on power.	Replace the battery in the charger. For instruc- tions on replacing your charger battery, see <i>In-</i> <i>stalling a battery in the charger, on page 5</i> .
During charging, the flashing green light on the charger turns off	Your transmitter is low on power.	 Charge the transmitter continuously for one hour. If flashing does not stop, pro- ceed to step 2.
and you see a series of quick flashing red lights on the charger for two seconds at a time.		2 Charge the transmitter continuously for eight hours. If flashing does not stop, call the 24 Hour HelpLine. It may be time to replace your transmitter.
During charging, a mix of quick and long flash- ing red lights appear on the charger.	Your charger and your transmitter are low on power.	1 Replace the battery in the charger. For in- structions on replacing your charger bat- tery, see <i>Installing a battery in the</i> <i>charger, on page 5.</i>
		2 Charge the transmitter continuously for one hour. If the quick flashing red lights do not stop, proceed to step 3.
		3 Charge the transmitter continuously for eight hours. If flashing does not stop, call the 24 Hour HelpLine. It may be time to replace your transmitter.

Problem	Likely Cause(s)	Resolution
The green light on the transmitter does not	Your transmitter is not fully connected.	1 Disconnect the transmitter from the sen- sor.
flash when you con- nect it to the sensor.	Your transmitter is low on power. Your sensor is not prop- erly inserted into your body.	2 Wait for five seconds and reconnect them. If the green light still does not flash, pro- ceed to step 3.
		3 Fully charge the transmitter and connect it to the tester. If the green light still does not flash, see troubleshooting on "The green light on the transmitter does not flash when you connect it to the tester". If the green light flashes, proceed to step 4.
		4 Disconnect the transmitter from the tester, wait at least five seconds, and connect the transmitter to the sensor. If the green light still does not flash, proceed to step 5.
		5 The sensor may not be properly inserted into your body. Remove the sensor from your body and insert a new sensor.
The green light on the transmitter does not flash when you con-	Your transmitter is not fully connected. Your transmitter is low on power.	1 Check the connection between the trans- mitter and the tester. If the green light still does not flash, proceed to step 2.
nect it to the tester.		2 Fully charge the transmitter.
		3 Test the transmitter with the tester again. If you still do not see the green light flash, call the 24 Hour HelpLine. It may be time to replace your transmitter.
Your transmitter bat- tery does not last for seven days.	tery does not last for fully charged when you	1 Fully charge the transmitter before con- necting it to the sensor. If the transmitter battery still does not last for seven days, proceed to step 2.
		2 Move away from any equipment that can cause RF interference. For more informa- tion on RF interference, see <i>Radio Fre-</i> <i>quency (RF) communication, on page 3</i> .
		3 Make sure your pump and your transmitter are located on the same side of your body to minimize any RF interference. If your fully charged transmitter battery continues to lose power before a full seven days, call the 24 Hour HelpLine. It may be time to replace your transmitter.

Problem	Likely Cause(s)	Resolution
Your transmitter has lost connection with your pump. Note: An alarm or alert occurs and a mes-	Your pump is out of range. There is RF interference from other devices.	1 Move away from any equipment that can cause RF interference. For more informa- tion on RF interference, see <i>Radio Fre-</i> <i>quency (RF) communication, on page</i> 3. If your transmitter is still not communicating with your pump, proceed to step 2.
sage appears when your transmitter has lost connection with your pump.		2 Make sure your pump and your transmitter are located on the same side of your body to minimize any RF interference. If your transmitter is still not communicating with your pump, call the 24 Hour HelpLine for assistance.

Storing the devices

Store the transmitter, charger, and tester in a clean, dry location at room temperature. If the transmitter is not in use, you must charge the transmitter at least once every 60 days. Although not required, you may store the transmitter on the charger. If you are storing the transmitter on the charger, you must disconnect and reconnect the charger and the transmitter at least once every 60 days. If the transmitter is left on the charger for more than 60 days, the battery will be permanently damaged.

Disposal

Discard the transmitter according to local regulations for battery disposal, or contact your healthcare professional for disposal information.

Biocompatibility	Transmitter: Complies with EN ISO 10993-1	
Applied parts	Transmitter	
	Sensor	
Operating conditions	Transmitter temperature: 32 °F to 113 °F (0 °C to 45 °C)	
	Caution: When operating the transmitter on a tester in air temperatures greater than 106 °F (41 °C), the temperature of the transmitter may exceed 109 °F (43 °C). This can cause the device to malfunction.	
	Transmitter relative humidity: 10% to 95% with no condensation	
	Transmitter pressure: 8.4 psi to 15.4 psi (57.6 kPa to 106 kPa)	
	Charger temperature: 50 °F to 104 °F (10 °C to 40 °C)	
	Charger relative humidity: 30% to 75% with no condensation	

Specifications

Storage conditions	Transmitter temperature: -4 °F to 131 °F (-20 °C to 55 °C)	
	Transmitter relative humidity: up to 95% with no condensation	
	Transmitter pressure: 8.4 psi to 15.4 psi (57.6 kPa to 106 kPa)	
	Charger temperature: 14 °F to 122 °F (-10 °C to 50 °C)	
	Charger relative humidity: 10% to 95% with no condensation	
Battery life	Transmitter: Seven days of continuous glucose monitoring immediately fol- lowing a full charge.	
	Charger: The charger uses one new AAA battery to charge the transmitter.	
Transmitter	2.4 GHz, 2M65G1D modulation, less than 1 mW ERP	
frequency		
Maximum output power (EIRP)	0.1 mW (-9.9 dBm)	
Radio Frequency (RF) communications	Pump to transmitter frequency: 2.4 Ghz; proprietary Medtronic protocol; range up to 6 feet (1.8 meters)	
	Utilizes the IEEE 802.15.4 protocol with proprietary data format.	
Transmitter expected service life	The transmitter expected service life is one year depending on patient us- age.	

Transmitter wireless communication

Quality of service

The transmitter and the MiniMed 670G insulin pump are associated as part of an 802.15.4 network for which the pump functions as the coordinator and the transmitter as an end node. In an adverse RF environment the pump will assess channel changing needs based on "noise" levels detected during an energy scan. The pump will perform the energy scan if after 10 minutes no CGM transmitter signal has been received. If the channel change occurs the pump will send beacons on the new channel.

The transmitter will initiate a channel search when beacon detection fails on the associated channel. The search will be conducted across all five channels. When the beacon is located the transmitter will rejoin on the identified channel. Upon re-association any missed packets (up to 10 hours) will be transmitted from the transmitter to the pump.

In normal operation the transmitter will transmit a packet every 5 minutes and retransmit the packet if the data is corrupted or missed.

Data security

The MiniMed 670G insulin pump is designed to only accept radio frequency (RF) communications from recognized and linked devices (you must program your pump to accept information from a specific device).

The MiniMed 670G insulin pump and system components (meters and transmitters) ensure data security via proprietary means and data integrity using error checking processes, such as cyclic redundancy checks.

Traveling by air

Your transmitter is safe for use on commercial airlines. If questioned by airline personnel about the use of your device, please show them your Medical emergency card. If they request that you turn off your system, you must comply.

Guidance and manufacturer's declaration

Guidance and Manufacturer's Declaration - Electromagnetic Emissions

The transmitter is intended for use in the electromagnetic environment specified below. The customer or the user of the transmitter should make sure that it is used in such an environment.

Emissions Test	Compliance	Electromagnetic Environment - Guidance
RF emissions CISPR 11	Group 1	The transmitter must emit electromagnetic energy in order to per- form its intended function. Nearby electronic equipment may be affected.
RF emissions CISPR 11	Class B	The transmitter is suitable for use in all establishments, including domestic and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.

Guidance and Manufacturer's Declaration - Electromagnetic Immunity

The transmitter is intended for use in the electromagnetic environment specified below. The customer or the user of the transmitter should assure that it is used in such an environment.

Immunity Test	IEC 60601 Test Level	Compliance Level	Electromagnetic Envi- ronment - Guidance
Electrostatic discharge (ESD) IEC 61000-4-2	±2 kV, ±4 kV, ±8 kV, ±15 kV Air ±2 kV, ±4 kV, ±6 kV, ±8 kV Contact	±2 kV, ±4 kV, ±8 kV, ±15 kV Air ±2 kV, ±4 kV, ±6 kV, ±8 kV Contact	For use in a typical do- mestic, commercial, or hospital environment.
Electrical fast transi- ent/burst	±2 kV for power supply lines	Not applicable	Requirement does not apply to this battery pow- ered device.

Guidance and Manufacturer's Declaration - Electromagnetic Immunity

The transmitter is intended for use in the electromagnetic environment specified below. The customer or the user of the transmitter should assure that it is used in such an environment.

Immunity Test	IEC 60601 Test Level	Compliance Level	Electromagnetic Envi- ronment - Guidance
IEC 61000-4-4	±1 kV for input/output lines		
Surge IEC 61000-4-5	±1 kV line(s) to line(s) ±2 kV line(s) to earth	Not applicable	Requirement does not apply to this battery pow- ered device.
Voltage dips, short in- terruptions and voltage variations on power supply lines	<5% <i>U</i> _T (>95% dip in <i>U</i> _T) for 0.5 cycle	Not applicable	Requirement does not apply to this battery pow- ered device.
IEC 61000-4-11			
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	3 A/m	400 A/m	Power frequency mag- netic fields should be at levels characteristic of a typical location in a typi- cal domestic, commer- cial, or hospital environment.

Note: U_T is the a.c. mains voltage prior to application of the test level.

Guidance and Manufacturer's Declaration - Electromagnetic Immunity

The transmitter is intended for use in the electromagnetic environment specified below. The customer or user of the transmitter should assure that it is used in such an electromagnetic environment.

Immunity Test	IEC 60601 Level	Compliance Lev- el	Electromagnetic Environ- ment Guidance
Conducted RF IEC 61000-4-6	3 V/m 150 kHz to 80 MHz	Not applicable	Not applicable

Guidance and Manufacturer's Declaration - Electromagnetic Immunity

The transmitter is intended for use in the electromagnetic environment specified below. The customer or user of the transmitter should assure that it is used in such an electromagnetic environment.

Immunity Test	IEC 60601 Level	Compliance Lev- el	Electromagnetic Environ- ment Guidance
Radiated RF IEC 61000-4-3	3 V/m 80 MHz to 2.5 GHz	10 V/m 80 MHz to 6 GHz	Portable and mobile RF com- munications equipment should be used no closer to any part of the transmitter, including ca- bles, than the recommended separation distance calculated from the equation applicable to the frequency of the transmit- ter.
			Refer to the recommended separation distance table for more information.
			d = 0.35 √P 80 MHz to 800 MHz
			d = 0.70 √P 800 MHz to 6 GHz
			Where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in meters (m).
			Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey ^a , should be less than the compliance level in each frequency range ^b .
			Interference may occur in the vicinity of equipment marked with the following symbol:

Guidance and Manufacturer's Declaration - Electromagnetic Immunity

The transmitter is intended for use in the electromagnetic environment specified below. The customer or user of the transmitter should assure that it is used in such an electromagnetic environment.

Immunity Test	IEC 60601 Level	Compliance Lev-	Electromagnetic Environ-
		el	ment Guidance

Note: At 80 MHz and 800 MHz, the higher frequency range applies.

Note: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption, and reflection from structures, objects and people.

^aField strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcasts and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the transmitter is used exceeds the applicable RF compliance level above, the transmitter should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating the transmitter.

^bOver the frequency range 150 kHz to 80 MHz, field strengths should be less than 3 V/m.

Recommended separation distances between portable and mobile RF communications equipment and the transmitter

The transmitter is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the transmitter users can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment and the transmitter as recommended below, according to the maximum output power of the communications equipment.

Rated maximum output power of transmitter (W)	Separation distance according to the frequency of transmitter (m)			
	150 kHz to 80 MHz Not applicable	80MHz to 800MHz d = 0.35 √P	800MHz to 6.0GHz d = 0.70 √P	
0.01	Not applicable	0.035	0.07	
0.1	Not applicable	0.11	0.22	
1	Not applicable	0.35	0.7	
10	Not applicable	1.1	2.2	
100	Not applicable	3.5	7	

For transmitters rated at a maximum output power not listed above, the recommended separation distance d in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where p is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.

English

Recommended separation distances between portable and mobile RF communications equipment and the transmitter

The transmitter is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the transmitter users can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment and the transmitter as recommended below, according to the maximum output power of the communications equipment.

Rated maximum output power of transmitter (W)	Separation distance according to the frequency of transmitter (m)			
	150 kHz to 80 MHz Not applicable	80MHz to 800MHz $d = 0.35 \sqrt{P}$	800MHz to 6.0GHz d = 0.70 √P	

Note: At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.

Note: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects, and people.

Warranty

Medtronic MiniMed, Inc. (or such other legal entity as may be referred to as manufacturer on the labeling of this device "Medtronic Minimed") warrants the Medtronic transmitter to the purchaser of the product against defects in material and workmanship for a period of one (1) year and the charger for up to one (1) year from the date of purchase.

During the warranty period, Medtronic MiniMed will replace or repair, at its discretion, any defective transmitter or charger, subject to the conditions and exclusions stated herein. This warranty applies only to new devices. In the event a transmitter or charger is replaced, the warranty period will not be extended past its original expiration date.

This warranty is valid only if the Medtronic transmitter or charger is used in accordance with the manufacturer's instructions. Without limitation, this warranty will not apply:

- If damage results from changes or modifications made to the transmitter or charger by the user, or third persons, after the date of purchase;
- If damage results from service or repairs performed by any person or entity other than the manufacturer;
- If damage results from a *Force Majeure* or other event beyond the control of the manufacturer;
- If damage results from negligence or improper use, including but not limited to: improper storage, submersion in water, physical abuse, (such as dropping);

• If damage results from use of the device in a manner other than according to the manufacturer's product labeling, instructions for use, or regulatory notifications.

This warranty shall be personal to the original purchaser. Any sale, rental or other transfer or use of the product covered by this warranty to or by a user other than the original purchaser shall cause this warranty to immediately terminate. This warranty does not apply to Glucose Sensors and other accessories.

The remedies provided for in this warranty are the exclusive remedies available for any breach hereof. Neither Medtronic MiniMed nor its suppliers or distributors shall be liable for any incidental, consequential, or special damage of any nature or kind caused by or arising out of a defect in the product.

All other conditions and warranties, other than mandatory statutory warranties, expressed or implied, are excluded, including the warranties of merchantability and fitness for a particular purpose.

This warranty gives the purchaser specific legal rights, and the purchaser may also have other rights that vary under local law. This warranty does not affect the purchaser's statutory rights.

Icon Table

SN	Serial number	
REF	Catalogue or model number	
(1x)	One transmitter, charger, and serter per container/package	
(2x)	Two testers per container/package	
M	Date of manufacture (YYYY-MM-DD)	
	Manufacturer	
8	Must refer to instruction manual before every use (appears blue on label).	
-20 °C -4 °F	Temperature limit	
((° i))	Non-ionizing electromagnetic radiation (RF communication).	
CONF	Configuration or unique version identifier	
Ŕ	Degree of protection against electric shock: Type BF applied part	

IP48	Transmitter: 4 is the level of protection against solid objects with a diameter above 1mm. 8 is the level of protection against the effects of continuous immersion in water [8 feet (2.4 meters) immersion for 30 minutes].
95%	Humidity limitation
R _{k Only}	Prescription only
Ţ	Fragile, handle with care
Ť	Keep dry
8	Recycle cardboard, paper, plastic packaging supplies and unwanted written material.
X	WEEE Initiative: DO NOT THROW IN TRASH. Recycle device according to local disposal requirements.
(NR)	Magnetic Resonance (MR) unsafe: keep away from magnetic resonance imaging (MRI) equipment.

©2016 Medtronic MiniMed, Inc. All rights reserved.

MiniMed[®] and Guardian[™] are trademarks of Medtronic MiniMed, Inc.

Bluetooth[®] is a registered trademark of Bluetooth SIG, Inc.

Ivory[®] is a registered trademark of The Procter & Gamble Company.

Detachol® is a registered trademark of Ferndale Laboratories Inc.

-26-





Medtronic MiniMed 18000 Devonshire Street Northridge, CA 91325 USA 800 646 4633 818 576 5555 www.medtronicdiabetes.com R Only MMT-7811 MMT-7736L MMT-7715 MMT-7512

6025977-013_a